

Cone Beam Volume Imaging (CBVI): 3D Applications for Dentistry

Robert A. Danforth, DDS

Diplomate, American Board of Oral and Maxillofacial Pathology
Associate Professor of Clinical Sciences, University of Nevada Las Vegas
School of Dental Medicine, Las Vegas NV.

Dale Miles, DDS, MS, FRCD(C)

Diplomate, American Board of Oral and Maxillofacial Radiology
Adjunct Professor, Arizona School of Dentistry & Oral Health and the University of Texas at San Antonio.

Impact upon Dental Imaging

The hallmark of dentistry in the new millennium is “**computer assisted dentistry**” and nothing characterizes this more than the impact 3D Cone Beam Computed Volume Imaging/Tomography (CBVI/CBVT) has had upon diagnosis, treatment planning and direct delivery of patient dental care. In less than a decade, CBVI has revolutionized oral and maxillofacial radiology. Although development was initially directed towards multiplanar viewing for dental implant and orthodontic treatment planning, secondary applications in other areas continue to expand because dentists find they are no longer confined to 2D, grayscale imaging and it provides information previously unobtainable (Table 1).

The intent of this article is to provide a brief overview of dental CBVI technology and examples of applications. While this technology has been almost universally described as CB Computed Tomography (CBCT), use of this term associated with the introduction of this technology created some concerns that this was in reality a medical device and, as such, should not be directly available to dentists for operation. This was particularly a concern in the USA and various terms such as CBVI or CB Volume Tomography (CBVT) were used to differentiate this connection with medical CT. In this article, CBVI is used to continue the differentiation.

Descriptions of CBCT devices for dentistry first appeared in the late 1990s with reports from Italy and Japan describing the development of prototype machines^{1,2}. The first devices to become commercially available were the Newtom QR 9000 (QR Italia, Verona, Italy) and the 3DX Accuitomo (J. Morita, Kyoto, Japan) which were initially introduced into the European and Asian markets. These were soon followed by the MercuRay (Hitachi, Japan, Twinsberg, OH, USA) and the i-CAT (Imaging Sciences, Hatfield, New Jersey, USA). All of the devices are capable of providing complete or partial skull fields of view (FOV) except for the

3DX Accuitomo which is a regional quadrant device with a field size of 40 x 30mm. Currently, there are approximately a dozen different machines either available, receiving government approval or in development (Table 2). Some of the recent machines have specifically targeted a limited FOV favoring regional, quadrant, and dental arch imaging similar to the initial Accuitomo device. Such an approach appears directed toward providing cost effective 3D imaging for the general practice market³⁻⁵.

How does it work? Why is it a better option for dentistry than medical Computed Tomography (CT)?

Most of the CBVI “stand-alone” devices which look like a panoramic machine and the patient either sits or stands during the examination. By contrast, the NewTom 9000, NewTom 3G, and the newer SkyView machines look like mini medical CT scanners and the patient lies in the supine position. However, the newest version of the Newtom VG marketed by Dent-X/AFP will be a stand-up version. While some companies have previously produced panoramic machines, currently, only the Planmeca ProMax CBVT is actually “upgradeable” from an existing panoramic platform

Depending upon the imaging device, image capture is generally a single 190⁰-360⁰ scan sweep around the patient’s head. The x-ray beam is “cone-beam” shaped aimed at a solid-state flat panel detector which covers the desired image volume in a single scan. The traditional “fan-shaped beam” employed in traditional medical CT requires multiple repeated slice scans to accomplish a similar image volume (Figures 1 A and B). This single scan CBVI approach is very successful at decreasing the patient’s absorbed x-ray dose from 6 to 15 times during the 3D imaging examination when compared to medical CT³.

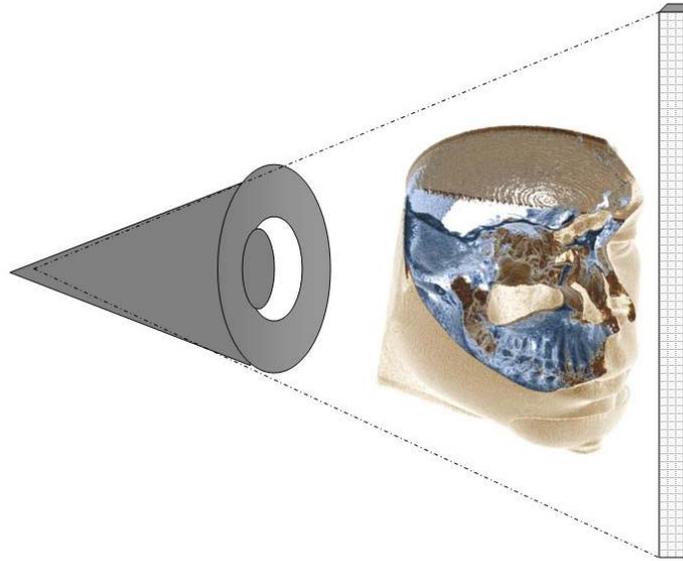


Figure 1 A

Cone-shaped beam rotates about head one time.



Figure 1 B

Fan-shaped beam rotates about patient's head hundreds of times acquiring thin horizontal or axial sections.

The cone beam scan results in several hundred dimensionally accurate isotropic images of the patient from different positions around the scan rotation. Data is immediately transferred to the computer which reconstructs it, in approximately 3-6 minutes depending upon the machine, into the anatomical volume for viewing at a precise 1:1 ratio with accurate dimensional measurements within 1/10 mm made across different tissue planes. Once reconstructed, multiplanar viewing of

the anatomical volume is accomplished with the imaging software. Volume data is saved in a DICOM (Digital Imaging and Communication in Medicine) file format to facilitate sharing images between clinicians and other imaging related services. Most manufacturers also provide a simpler “user” PC viewing program which contains basic 3D imaging tools. This allows referral sources to participate in limited 3 D imaging without having to purchase expensive software. For clinicians who wish to more actively plan their implant cases, third party software is available at a wide range of cost.

Applications and a new perspective for interpretation, diagnosis and treatment planning

The results of computer generated 3D visualization of multiplanar images are new applications which improve interpretation, diagnosis and treatment planning. Imaging for orthodontic and dental implants are prime examples. Full volume CBVI data for orthodontics has spawned new research in the orthodontic 1:1 ratio cephalometric measurements for treatment planning and several software programs have been developed for this purpose (Figure 2).

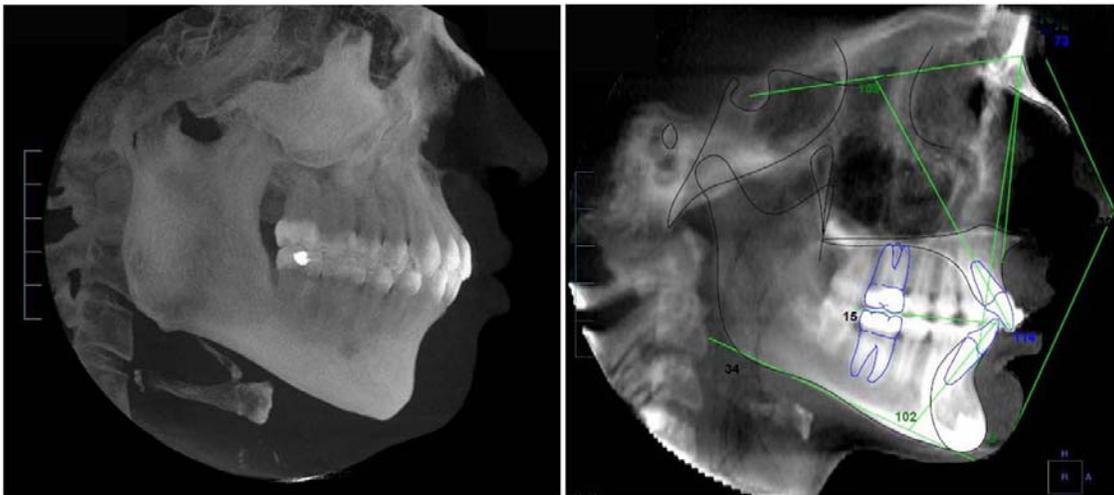


Figure 2A

Figure 2B

Lateral cephalometric images made using Multiple Image Planes (MIP). Figure 2A MIP images with soft tissue profile and Figure 2B image with orthodontic measurements using software developed for CBVI.

There has also been a tremendous impact upon dental implant treatment and treatment planning. Arch and ridge configuration, location of the nerve canals and sinus cavities can now be specifically identified prior to surgery thus decreasing surgical uncertainty and potential patient morbidity (Figures 3A and 3B).



Figure 3A

Analytical image to render inferior alveolar canal structures visible in colors and precise implant site measurements. The proposed implant fixture may now be placed without fear of impinging on the inferior alveolar nerve.

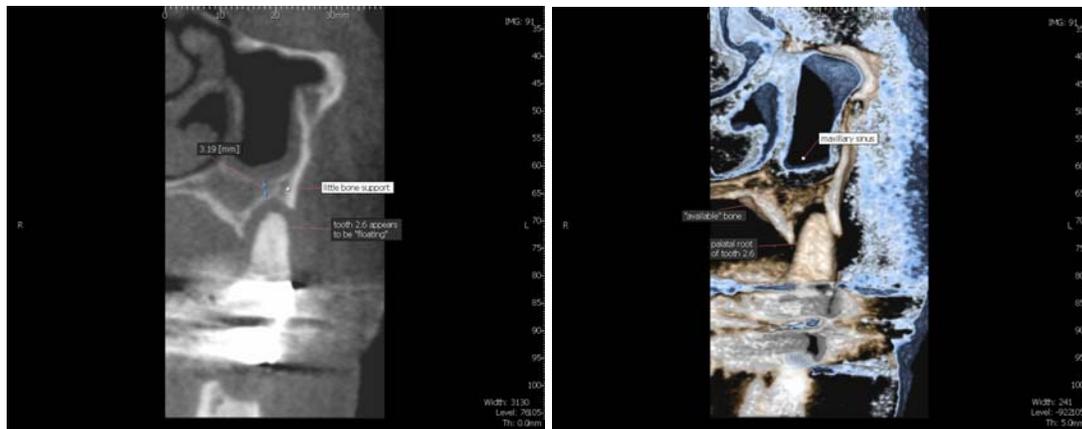


Figure 3B

Proposed implant site of tooth to be extracted showing insufficient bone, periapical disease and associated inflammatory changes in maxillary sinus. Clearly, there is other treatment necessary before even attempting placement of an endosseous implant in this site.

Besides mere visualization of cross-sectional imaging of jaw site locations, innovative software programs have been developed to make custom precision surgical guides, virtual models, and laser generated resin models, all of which advance treatment planning and assist delivery capabilities (Figure 4). And, the advantage of these innovative changes is that they are relatively cost effective and readily available to the dentist. Direct fabrication of dental appliances from CBVI data using computer assisted design-computer assisted manufacturing (CAD-CAM) systems would seem a next progressive development.

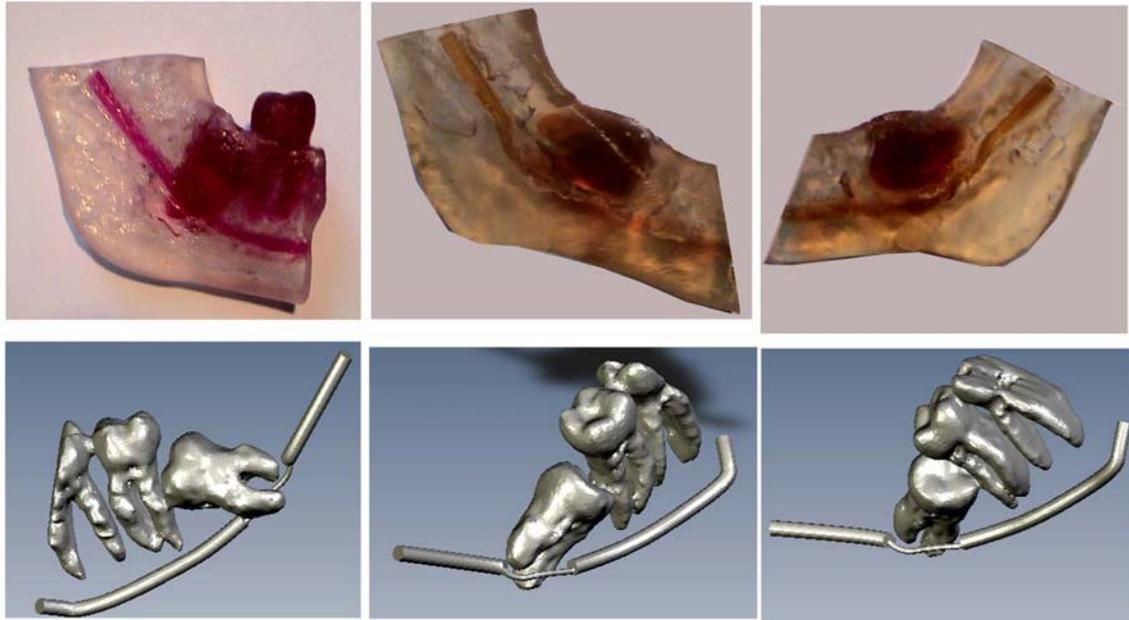


Figure 4.

Top row, reconstructed laser-generated resin models from CBVI volume images to assist surgical treatment planning for impacted third molars at risk for injury to the inferior alveolar nerve. Bottom row, visual models segmented from CBVI image data. Interactive tools allow for rotation, measurements, and image “slicing and dicing” to assist surgical planning.

Beyond dental- 3D image volume produces new responsibility,

Although CBVI has created exciting new 3D methods for dental applications, it has also created increased responsibility for image interpretation. Unlike superimposed 2-D imaging, a large volume of viewable collected patient data exists outside of the specific region of interest (ROI) in the multiplanar sections. Such image data must be viewed in the X, Y, and Z planes to ensure that a comprehensive evaluation has been performed and documented⁶. Failure to recognize these or to act upon them could result in various pathologic conditions going undetected that could have a negative outcome on either the provided treatment, or worse, patient health and longevity.

Table 1, part 2 indicates that sinus and airway are applications for evaluating the presence of conditions involving pathologic lesions, sleep apnea or head and neck artery calcifications. In some cases, these may be secondary/incidental but life-threatening findings which are most likely located outside the jaws in non-ROI 3D multi-planes. As such, many dentists may not be as familiar or comfortable with looking for or interpreting these conditions. Despite any unfamiliarity, many recent reports indicate dentistry has a role in such conditions as; airway sleep apnea problems, sinus conditions, bisphosphonate jaw osteonecrosis, and awareness of carotid and other arterial calcifications as an indicator of potential

stroke or advanced diabetic/renal problems ⁷⁻¹¹. These are examples of expanded regions of concern and responsibility where CBVI can assist detection of these conditions. In these circumstances, referral to a diagnostic imaging radiologist/specialist for an interpretative report of the 3D data volume is recommended (Figure 5 A – C).

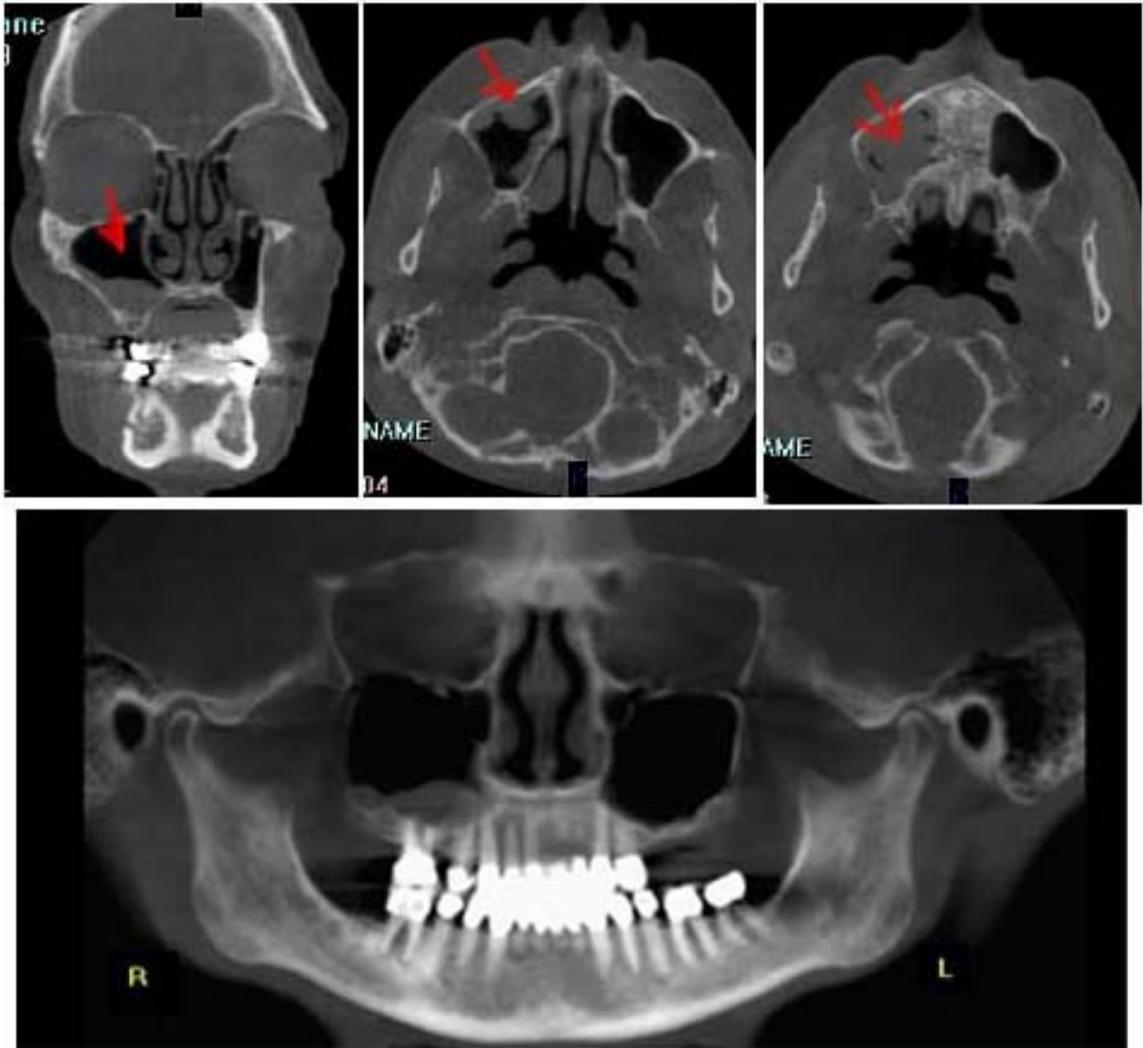


Figure 5 A

The extent of right maxillary sinus disease is shown using multiple plane reviewing and a reconstructed panoramic image. Use of such techniques better defines the location of the disease and any impact it may have on the proposed treatment plan.

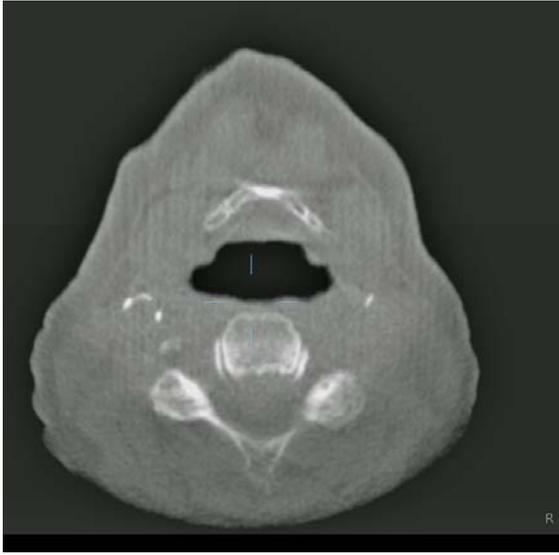


Figure 5B

Bilateral circular calcifications in the region of the carotid sheath at the level of C3/C4 consistent with MAC (medial arterial calcinosis) seen in diabetic patients especially with end-stage renal disease (ESRD).

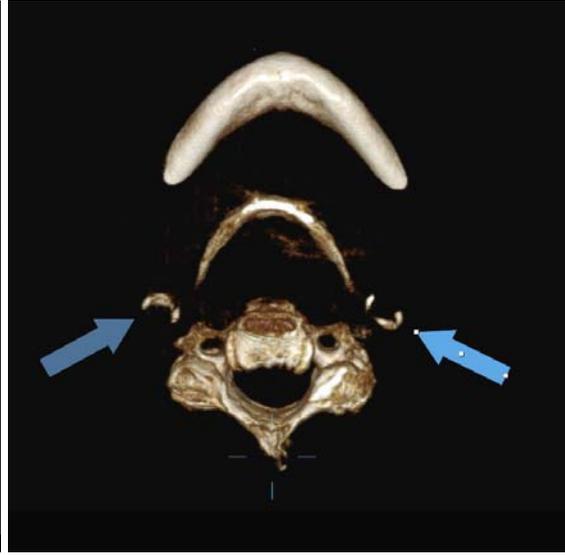


Figure 5C

Conclusion

Any clinician reviewing the Table 1 application list, knows the limitations imposed by 2D imaging in a 3D world when interpretive diagnostic dilemmas arise in these areas. The length of the list indicates that, actually, no aspect of dentistry is free from such dilemmas. Until CBVI, when confronted with such situations, clinicians either chose medical related CT imaging or continued to rely upon their tradition interpretation of 2D imaging despite inherent disadvantages. CBVI imaging provides relief for clinicians from 3D mental reconstructions by using outstanding computer assisted 3D imaging technology for multiplanar image visualization, diagnosis and treatment planning.

The CBVI data obtained is similar to that from a high-end medical CT scan, but at a lower financial cost to the patient and at much reduced radiation dose to the patient¹². Since these devices are becoming widely available in either dental offices or dental imaging centers, reliance upon experience based guess-estimations is being replaced accurate imaging technology. As such treatment outcomes are improved and both patients and dentistry benefit. Despite these obvious advantages, comprehensive 3D image evaluation associated with CBVI is a new responsibility for dentistry that must be accomplished in order to comply with the ethics of comprehensive patient care.

For more information on CBCT and volumetric images go to www.learn.digital.net

References:

1. Mozzo P, Proccacci A et al : A new volumetric CT machine for dental imaging based on the cone-beam technique: preliminary results. *Eur Radiol* 8:1558-64,1998)
2. Arai Y, Tammsalo E, et al: Development of a compact computed tomographic apparatus for dental use. *Dentomaxillofac Radiol* 28:245-8,1999)
3. Farman AG, Levato CM, and Scarfe WC. A primer on Cone Beam CT. Inside Dentistry, Vol. 1, 90-93,2007.
4. Miles DA and Danforth RA: A Clinician's Guide to Understanding Cone Beam Imaging, PennWell, (in press).
5. Hayakawa Y; Newly added oral and maxillofacial radiology imaging technologies. *AADMRT Newsletter*, winter 2007.
6. Miles DA: Clinical Experience with Cone Beam Volumetric Imaging: Report of Findings in 381 Cases, *US Dentistry*, 2006;1(1):39-41.
7. Friedlander AH, Garrett NR, Norman DC: The prevalence of calcified carotid artery atheromas on the panoramic radiographs of patients with type 2 diabetes mellitus. *JADA*, 2002;133:1516-23.
8. Miles DA and Craig RM: The calcified facial artery: A report of its panoramic appearance and incidence. *Oral Surg*1983; 55(2):214.
9. Almog DM, and et al: Correlating carotid artery stenosis detected by panoramic radiography with clinical relevant carotid artery stenosis determined by duplex ultrasound. *Oral Surg Oral Med Oral Path Oral Radiol Endod* 2002;94:768
10. Magliocca KR and Helman JI: Obstructive sleep apnea: Diagnosis, medical management and dental implication. *J AM Dent Assco*. 2005;136:1121-1129.
11. Kumar V, Pass B, et al: Bisphosphonate-related osteonecrosis of the jaws. *JADA*, 2007(138):602-609.
12. Ludlow JB, Davies-Ludlow LE, Brooks SL and Howerton WB. Dosimetry of 3 CBCT devices for oral and maxillofacial radiology: CB Mercuray, NewTom 3G, and i-CAT. *Dentomaxillofac Radiol* 2006 35:219-226.

Table 1. Current dental applications of CBVI

1. Primary applications

- Pre-surgical implant site assessment
 - Inferior alveolar nerve location
 - Alveolar ridge configuration
 - Spatial relationships of anatomical structures
 - Surgical guide construction

- Orthodontic assessment
 - Cephalometric and space 1:1 analysis

- Temporomandibular joint evaluation
 - Multiplanar views, Dentition relationship

- Extraction cases
 - Impactions, third molars/nerve location, canines
 - Supernumerary teeth

2. Secondary applications

- Paranasal sinus evaluation
- Airway studies for sleep apnea
- Visualization of jaw pathology
- Trauma evaluation
- Endodontic evaluation
- 3D virtual models
- TADs (Temporary Anchorage Devices)
- Other CAD/CAM devices

Table 2 CBVI devices and manufacturer: 2007

1. **NewTom 3G, VG:** Dent-X, Elmsford, New York, USA
2. **Accuitomo 3DX:** J. Morita, Kyoto, Japan
3. **i-Cat:** Imaging Sciences (Danaher owned), Hatfield, Pennsylvania, USA
4. **CB MercuRay, CB Throne:** Hitachi Medical Systems, Twinsberg, OH, USA
5. **Iluma,** Imtec Corporation (Kodak distribution), Ardmore, Oklahoma, USA
6. **Galileos:** Sirona Dental Systems, Bensheim, Germany
7. **ProMax 3D CBVT,** Planmeca, Helsinki, Finland
8. **SkyView:** MyRay Dental Imaging, Imola, Italy
9. **TeraRecon:** Yoshida Dental Manufacturing, Tokyo, Japan
10. **Alphard and PSR900N:**Asahi Roentgen, Kyoto, Japan
11. **Scanora 3D:** Soredex (PaloDex), Tuusula, Finland
12. **PicassoPro:** eWoo, (Gyeonggi-do, Korea)

Table 3 Viewing Software for Clinicians

SimPlant: Materialise, Leuven, Belgium

VIP (Virtual Implant Placement): Implant Logic, Systems, Cedarhurst, NY, USA

V Implant, CB Works, Accurex, On Demand 3D: CyberMed International, Seoul, Korea

In Vivo, Anatomage, San Jose, CA, USA

NobelGuide™: Nobel Biocare, Göteborg, Sweden